

Address

#204 B, Falcon Gallery Building,
Al Mina Road, Near Jumaira 1,
Bur Dubai, Dubai, UAE.
P.O.Box: 28842.

Phone

+971 4 354 2905

Email:

info@spectronixglobal.com

www.spectronixmedical.com



VITAL INJECTOR™

VITAL INJECTOR™




Enhance your results for more
effective drug delivery, PRP and
Mesotherapy

How does it work?

Vital Injector is a non-invasive intradermal injecting device using a syringe. It maximizes the treatment effect and shortens the procedure time by adopting 5 multi-Needles.

Technology

Vital Injector injects drug accurately and equally by adopting 5 multi-Needles. Its safe vacuum technique minimizes drug loss and alleviates pain by improving the speed of needling. It is the system with the convenient needles and a removable syringe.

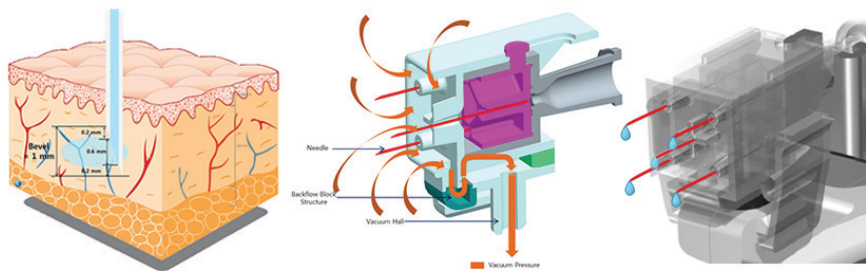

Components



Main body



Injector



Multi-needle(Option)

Specifications

Needle Information: 31 G, 5Pin Multi Needle Type

Needle Penetration Control: 0 ~ 5 mm (Interval 0.2mm)

Syringe Type : 1cc, 2cc, 3cc, 5cc

Operation Mode: Intermittent, Continuous

Voltage & Frequency: DC 12V or DC 24V, 50~60 Hz

Dimension (W x D x H mm) : Injector : 303 x 51 x 175

Body : 232.6 x 247.9 x 126 | **Panel :** Gun LCD : 1.4 Inch

Body LCD : 4.3 Inch, Touch Screen | **Weight :** 2.4 Kg (Body + Injector)



A . Design of Trigger

The trigger is designed to maintain the center of gravity of the device and inject steadily while handling for the treatment.

B . 5Multi-Needle

Maximizing the treatment effect by adopting 5 Multi-Needles! Convenient removable Multi- needles system, alleviating pain by improving the speed of needling.

C . Various types of Syringes are applicable

4 syringe types(1cc, 2cc, 3cc & 5cc) can be applied and highly compatible.

D . Display Panel

It shows Information for the treatment like the number of injection, the depth of needle penetration, drug dosage, and the operation mode, etc.

E . Move Forward & Back Button . Standby/Ready Button

Control button.

F . Position Holder

Stable performance of the syringe pusher by installing the screw shaped fastening device.

G . System cable

You can charge it on the provided cable mounting charger.

H .Main Body

The LCD touch screen enables you to control the number of injection, the depth of needle penetration, drug dosage, the operation mode, and the speed/vacuum Level.

I . Filter

Hygienic maintenance with regular changes of filters located in the lower part of the product.

VITAL INJECTOR™